



EN



New Frontier of Lithium Disilicate-Based
CAD/CAM Blocks & Disks

Amber[®] Mill

User's Manual



www.hassbio.com



RX Only

Human-Aid
System Supplier

beLIVE
HASS

Amber[®] Mill

사용자 메뉴얼

Table of Contents

1	Introduction	3
2	Preparation Guide	4
3	Block Selection	5
4	CAD/CAM Process	6
5	Before Crystallization	7
6	Crystallization for HT(High Translucency)	8
7	Crystallization for MT(Medium Translucency)	9
8	Crystallization for LT(Low Translucency)	10
9	Crystallization for MO(Medium Opacity)	11
10	Crystallization Results	12
11	Re-Firing	13
12	Co-Firing	14
13	Rapid Mode	15
14	Characterizing	16
15	Staining & Glazing	17
16	Cementation	18
17	Indications / Contra-Indications	19
18	Product Line-up	20
19	FAQs	21

1. Introduction



New Frontier of Lithium Disilicate-Based
CAD/CAM Blocks & Disks

Amber[®] Mill

Aesthetic values

Resulting from its excellent opalescence, fluorescence and innate controllable value, Amber[®] Mill restorations feature the most natural and lifelike looks in multi-chromatic gradations.


Structural stability

Cross-linked lithium disilicate after translucency heat-treatment allows that Amber[®] Mill crowns can achieve a high mechanical properties. (over 450 MPa)


Edge stability on NLD technology

Outstanding machinability of Amber[®] Mill is clearly proven by stable edges of the milled restorations. Less chipping occurrence with Amber[®] Mill tells you that Amber[®] Mill is the machinable lithium disilicate block for CAD/CAM system.


2. Preparation Guide

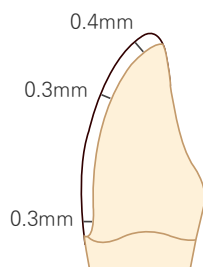


Thin Veneer

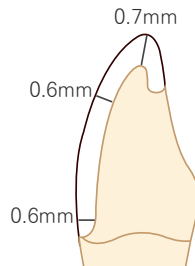


Veneer

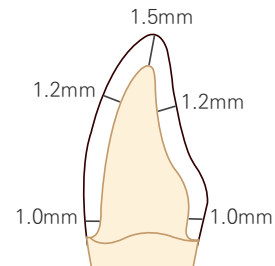

Anterior Crown

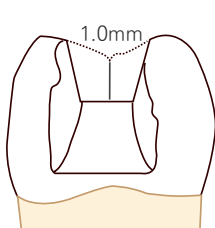

Inlays

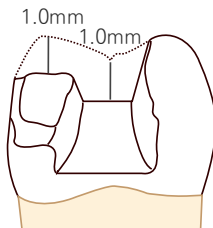

Onlays


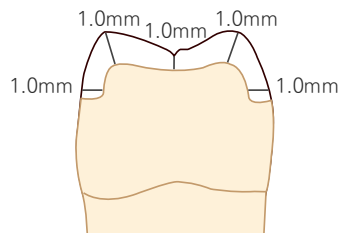



Table Top

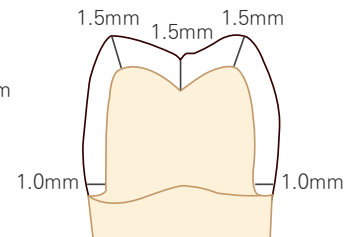


Posterior crown

3. Block Selection


Processing Technique & Indications

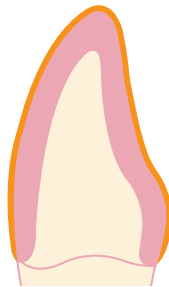

Staining technique
HT / MT / LT

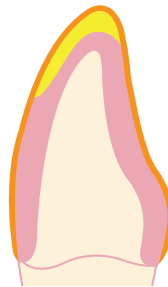

Cut-back technique
HT / MT / LT / MO

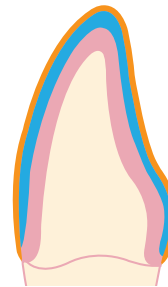

Layering technique
HT / MT / LT / MO

- Amber® Mill**
- staining & glazing
- incisal veneering material
- veneering material


Indications								
Table Tops	Thin Veneers	Veneers	Inlays	Onlays	Partial Crowns	Anterior Crowns	Posterior Crowns	3-Unit Bridges
High Translucency								
Medium Translucency								
Low Translucency								
Medium Opacity								

TIP!



For use in anterior and premolar bridge restorations can be expected.

4. CAD/CAM Process

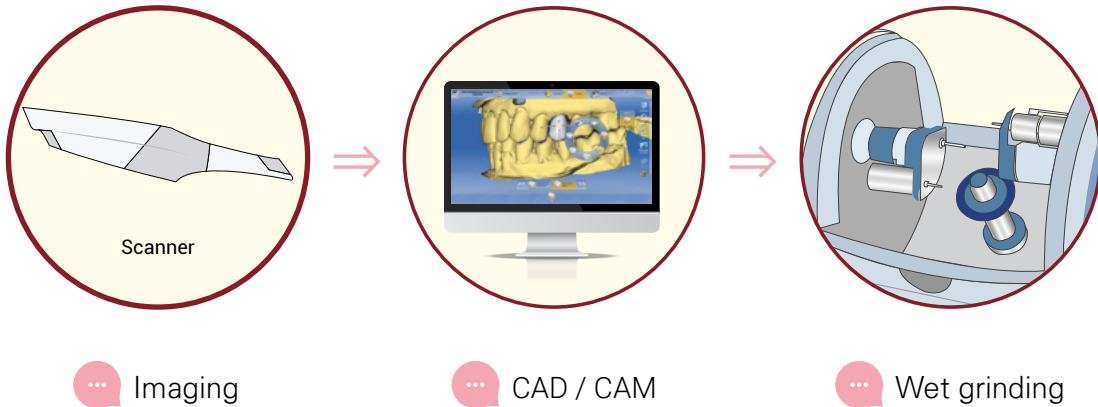

TIP!



Amber Mill is available in various CAD / CAM systems.




5. Before Crystallization



... Finish surface with grinding instrument.



... Use peg putty before crystallization.

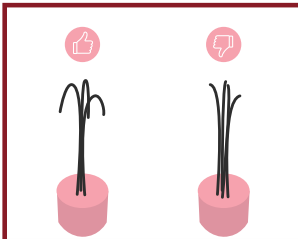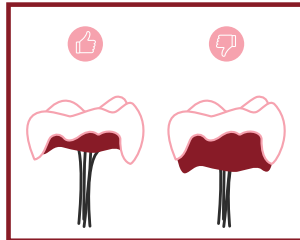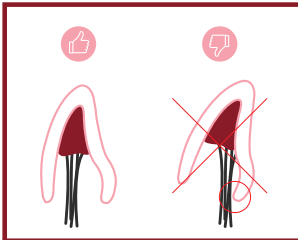

TIP!



Use peg putty on thin-metal pin as little as possible to minimize the absorption of heat by peg putty.

6. Crystallization for HT(High Translucency)

... DEKEMA Austromat 624i¹⁾

	Standard Mode			Rapid Mode		
Dry			--			--
Close			02:00			01:00
Preheat	450°C		01:00	450°C		01:00
Temperature 1	830°C	60°C /min	15:00	790°C	100°C /min	--
Temperature 2	690°C	60°C /min	--	830°C	15°C /min	05:00
Temperature 3	—°C	—°C /min	--	680°C	70°C /min	--
VAC (off/level/hold)	830°C	100%	15:00	830°C	80%	05:00



TIP!

Watch more detail video!

... IVOCLAR VIVADENT PROGRAMAT CS²⁾

Standard Mode

B °C	S min.	t / °C/min.	T °C		H min.	VAC. 1/ VAC. 2 °C		L °C	tL*
400	3.00	60	HT	815	15.00	HT	550/815	690	0

Rapid Mode

B °C	S min.	t ₁ ↙ °C/min.	t ₁ °C	H min.	t ₂ ↙ °C/min.	t ₂ °C		H min.	VAC. 1/ VAC. 2 °C			L °C	tL*
400	1.00	90	780	0:00	30	HT	815	3.00	HT	690/780	780/815	690	40

1) Austromat 624i is a registered trademark of DEKEMA.

2) PROGRAMAT CS is a registered trademark of IVOCLAR VIVADENT.

7. Crystallization for MT(Medium Translucency)

... DEKEMA Austromat 624i¹⁾

	Standard Mode			Rapid Mode		
Dry			--			--
Close			02:00			01:00
Preheat	450°C		01:00	450°C		01:00
Temperature 1	840°C	60°C /min	15:00	800°C	100°C /min	--
Temperature 2	690°C	60°C /min	--	840°C	20°C /min	05:00
Temperature 3	—°C	—°C /min	--	680°C	70°C /min	--
VAC (off/level/hold)	840°C	100%	15:00	840°C	80%	05:00



TIP!

!

Watch more detail video!

... IVOCLAR VIVADENT PROGRAMAT CS²⁾

Standard Mode

B °C	S min.	t / °C / min.	T °C		H min.	VAC. 1/ VAC. 2 °C		L °C	tL *
400	3.00	60	MT	825	15.00	MT	550/825	690	0

Rapid Mode

B °C	S min.	t ₁ ↙ °C/min.	t ₁ °C	H min.	t ₂ ↗ °C/min.	t ₂ °C		H min.	VAC. 1/ VAC. 2 °C			L °C	tL*
400	1.00	90	780	0:00	30	MT	830	3.00	MT	690/780	780/830	690	40

1) Austromat 624i is a registered trademark of DEKEMA.
2) PROGRAMAT CS is a registered trademark of IVOCLAR VIVADENT.

8. Crystallization for LT(Low Translucency)

... DEKEMA Austromat 624i¹⁾

	Standard Mode			Rapid Mode		
Dry			--			--
Close			02:00			01:00
Preheat	450°C		01:00	450°C		01:00
Temperature 1	855°C	60°C /min	15:00	800°C	100°C /min	--
Temperature 2	690°C	60°C /min	--	855°C	20°C /min	04:00
Temperature 3	—°C	—°C /min	--	680°C	70°C /min	--
VAC (off/level/hold)	855°C	100%	15:00	855°C	80%	04:00



TIP!
! Watch more detail video!

... IVOCLAR VIVADENT PROGRAMAT CS²⁾

Standard Mode

B °C	S min.	t / °C / min.	T °C		H min.	VAC. 1/ VAC. 2 °C		L °C	tL*
400	3.00	60	LT	840	15.00	LT	550/840	690	0

Rapid Mode

B °C	S min.	t ₁ ↙ °C/min.	t ₁ °C	H min.	t ₂ ↗ °C/min.	t ₂ °C		H min.	VAC. 1/ VAC. 2 °C			L °C	tL*
400	1.00	90	780	0:00	30	LT	845	3.00	LT	690/780	780/845	690	40

1) Austromat 624i is a registered trademark of DEKEMA.

2) PROGRAMAT CS is a registered trademark of IVOCLAR VIVADENT.

9. Crystallization for MO(Medium Opacity)

... DEKEMA Austromat 624i¹⁾

	Standard Mode			Rapid Mode		
Dry			--			--
Close			02:00			01:00
Preheat	450°C		01:00	450°C		01:00
Temperature 1	875°C	60°C /min	15:00	800°C	100°C /min	--
Temperature 2	690°C	60°C /min	--	870°C	25°C /min	04:00
Temperature 3	—°C	—°C /min	--	680°C	70°C /min	--
VAC (off/level/hold)	875°C	100%	15:00	870°C	80%	04:00



TIP!



Watch more detail video!

... IVOCLAR VIVADENT PROGRAMAT CS²⁾

Standard Mode

B °C	S min.	t °C/min.	T °C		H min.	VAC. 1/ VAC. 2 °C		L °C	tL*
400	3.00	60	MO	860	15.00	MO	550/860	690	0

Rapid Mode

B °C	S min.	t ₁ ↙ °C/min.	t ₁ °C	H min.	t ₂ ↙ °C/min.	t ₂ °C		H min.	VAC. 1/ VAC. 2 °C			L °C	tL*
400	1.00	90	780	0:00	30	MO	865	3.00	MO	690/780	780/865	690	40

1) Austromat 624i is a registered trademark of DEKEMA.


2) PROGRAMAT CS is a registered trademark of IVOCLAR VIVADENT.

10. Crystallization Results



Temperature

Low Temp. —————> High Temp.

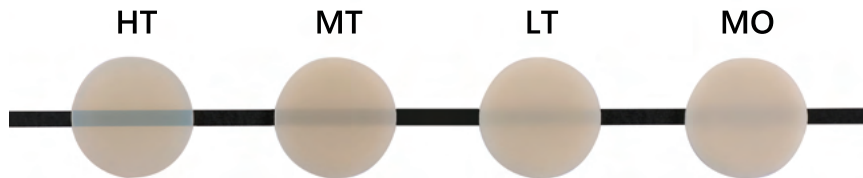

Translucency




Four different translucencies (HT, MT, LT, MO) in one block.

Choose your own translucency just by changing crystallization temperature.

11. Re-Firing

... Change of Opacity(translucency) by Re-Firing

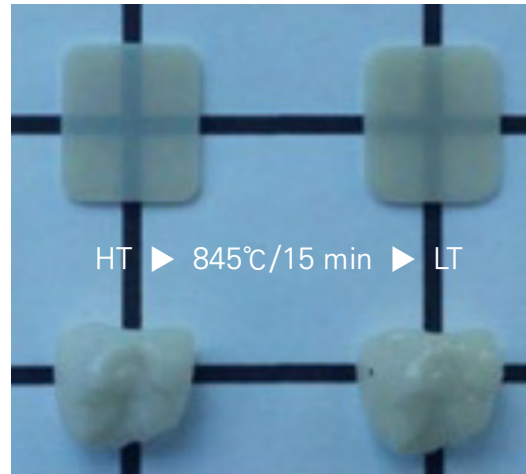


Amber® Mill A2

TIP!




For value color, it is achievable to lower their translucency by re-firing them. For example, you may apply **5°C higher** re-firing than normal low translucency(LT) temperature to high translucency(HT) crowns and keep the same holding time of 15 minutes so that the final crowns can be low translucency(LT).

12. Co-Firing

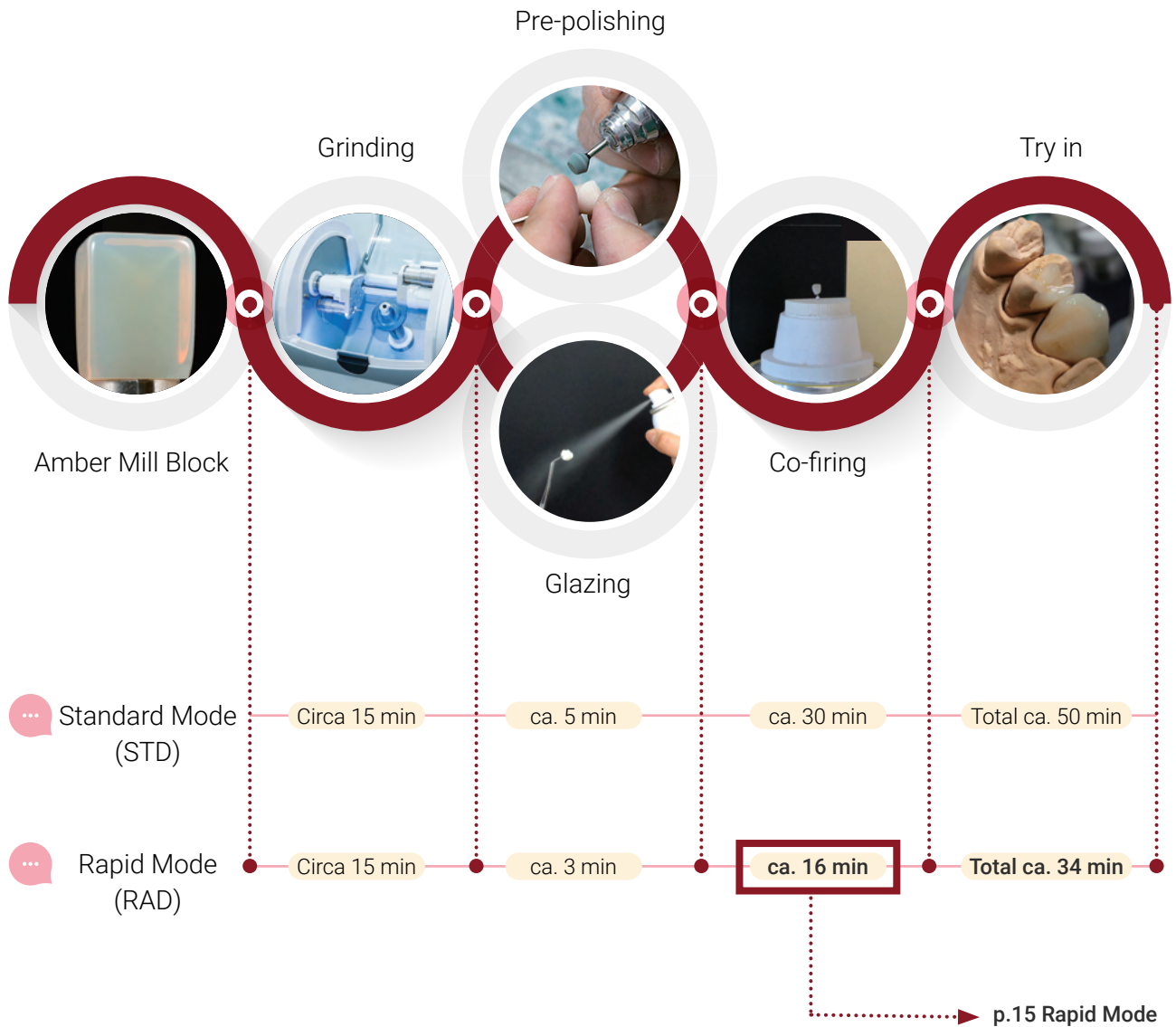


13. Rapid Mode

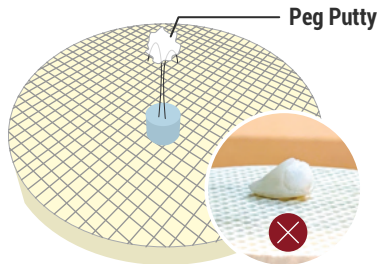



How to apply peg putty

Instead of filling up the inside of the restoration with peg putty, we recommend you to put the minimum amount of peg putty on the metal pin and attach restoration to the pin.

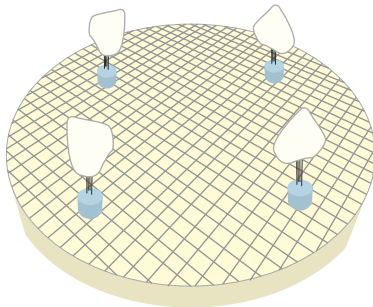


Quantity and arrangement method

When arranging multiple restorations, equal spacing should be used between all restorations.
(Ex. 4 restorations, it is ideal to place them in 4 equal directions)

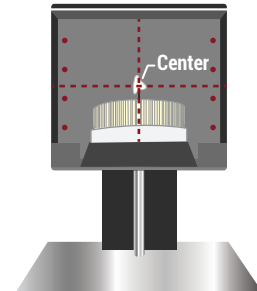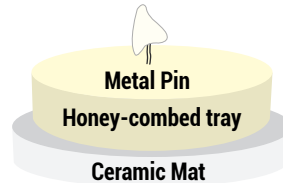

Placement of restorations in the furnace

Adjust the height of the pin to be in the center of oven chamber. This will assure constant and equal heat transfer to the restoration.

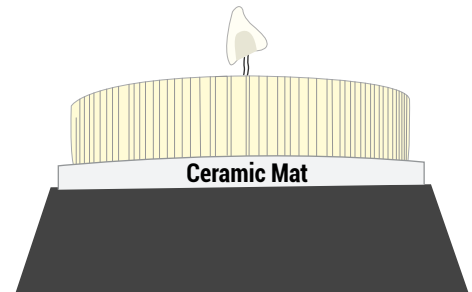

Placement of ceramic mat

A Ceramic mat is required to be placed under honey-combed tray to minimize heat loss.

RAPID MODE



TIP!

! Recommendation : **Single Unit** and/or **Spray type glaze**.

! This is based on Austromat 624i of DEKEMA.

14. Characterizing



... Ensure the surface is clean by removing bubbles. At this time, it is also necessary to work while applying water.

15. Staining & Glazing



- ... Apply the stain in accordance with the target shade. Select Stain & Glaze products that can be fired at temperature lower than 860°C.

16. Cementation



... After fit-in, etch with Hydrofluoric acid 4% for 30 seconds.



... Rinse with water and dry with an air syringe.



... Silanize the inner aspects for 20 seconds and dry with an air syringe.



... Adhesively lutes the prosthesis with self-adhesive resin cement.

TIP!



Please keep the etching time indication; exceeding the time can cause fragility.



Please refer to the manufacturer's guideline for the usage of silane.

17. Indications / Contra-Indications

Indications



Inlays



Onlays



Veneers



Anterior Single Crowns



Posterior Single Crowns



3-Unit Bridge
*up to the second Premolar

Contraindications

- Very deep subgingival preparations
- Maryland bridges
- Patients with severely reduced residual dentition
- Bruxism
- Cantilever bridges

19. FAQs

Q 1. What is Amber[®] Mill's strength?(Before & after the translucency heat-treatment)



Q 2. What powders are compatible with Amber[®] Mill?

A Amber[®] Mill is compatible with a wide variety of veneering powders. As to the powders for lithium disilicate, those powders with CTE(coefficient of thermal expansion) less than or equal to $10.0 \times 10^{-6}/^{\circ}\text{C}$ are compatible. Zirconia powders with baking temperature under 850°C are also compatible with Amber[®] Mill.

Q 3. How is the translucency heat-treatment temperature difference changing the translucency?

A When applying the translucency heat-treatment, distribution of fine crystalline and coarse crystalline as well as crystal density are changed, which changes the transmission ratio of visible light. Usually, the higher heat-treatment temperature gets, the more coarse crystalline and density Amber[®] Mill block has. This scatters light ray more and, consequently, the translucency gets lower.

Q 4. What should be mainly considered for the translucency heat-treatment?

A Combination of two factors-temperature and holding time-for the translucency heat-treatment of Amber[®] Mill differentiates the resulted translucency. Based on the recommended translucency heat-treatment schedule, each user is advised to verify his or her own optimized conditions for the furnace to use. Once the optimized version is identified, you will be able to create a wide range of translucency with just one Amber[®] Mill block and choose the exact translucency level as targeted.

Q 5. Is it possible to change the translucency by re-firing?

A For highly translucent restorations, it is achievable to lower their translucency by re-firing them. For example, you may apply 5°C higher heating than normal low translucency (LT) translucency temperature to high translucency (HT) crowns and keep the same holding time of 15 minutes so that the final crowns can be low translucent (LT).

Q 6. What about the chemical durability?

A When ceramic materials are exposed to moisture for a long time, degradation of material properties usually takes place because of ionic dissociation. Production process of Amber[®] Mill hires high purity raw materials and minimizes uninvited alkali impurities so that it generates fine crystalline structure and chemically durable glass matrix. This is why Amber[®] Mill performs more excellent chemical durability than existing products and ensures high long-term clinical reliability.



HASS Corporation

77-14, Gwahakdanji-ro, Gangneung-si, Gangwon-do, KOREA 25452

Tel: +82-70-7712-1300 / Fax: +82-33-644-1231

Customer Support : +82-2-2083-1367

E-mail : hasscorp@hassbio.com

Website : www.hassbio.com

This material is designed for usage in dentistry. Follow instructions HASS is not liable for any loss caused by failure to comply with regulation or scope of indication.

Users are responsible for testing products to verify the compatibility for any usage which are not written in the instructions. The explanations and data contained within do not carry any guarantees and/or obligations.

All enclosed recommendations and restrictions apply when used with products from other manufacturers.